



MOUNTING REFERENCE GUIDE

VESA 75/100 Adaptation

CIM Cable Integrated Mounting System



THINK BEFORE
YOU PRINT



Content

General notes	3
How to create part numbers	4
Product highlights	6
Flush mount	7
Support arm pivot	9
Extended pivot arm light duty	11
Extended pivot arm heavy duty	13
Height adjustable arm	15
Extended height adjustable arm	17
High position height adjustable arm	19
Extended high position height adjustable arm	21
Counter top	23
Counter top for deep monitors	25
Cart light duty	27
Cart light duty kit	29
Fixings	31

General notes

Monitor adaptations

VESA adaptations are available in various versions:

- Standard VESA 75/100
- VESA rotatable (portrait/ landscape applications)
- VESA quick release for easy disassembly - refer to catalogue VESA 75/100 with quick release

Components' grounding

All arms with integrated cable management are equipped with **component grounding**. No grounding is required for arms with external or semi-integrated cable management with cables not being led through joints.

Load capacity

All **height adjustable arms** have a load capacity of 22 kg / 48.5 lbs. Rear mounted monitor adaptations and horizontally mounted 5" plunge plates accept a maximum monitor weight of 18 kg / 39.7 lbs. An additional 4 kg / 8.8 lbs may be mounted on a down post underneath the front swivel part. Depending on material and tensile factor all weight indications have a 4-fold and/ or 6-fold safety factor according to IEC 60601-1 standards.

Custom-designed solutions and MRI applications

Please contact us for **custom-designed** solutions.

Product marking

All products which are classified as medical product risk class 1 are CE marked. As **medical product risk class 1** they are conform with the regulations (EU) 2017/745 (MDR) and IEC 60601-1. Other standards of the IEC 60601 series are not applicable to our products, regardless of their version.

MD A stand-alone product that is already classified as a medical product class 1 before being integrated into an overall medical device. The product is CE-marked. A declaration of conformity is available.

NO Accessories or spare parts that have no independent function and are therefore not subject to the MDR. These products are not CE-marked. A declaration of conformity is not required.

Warranty

All CIM products have a 5-year warranty to be free of defects in material and functionality from the date of delivery.



How to create part numbers

Part numbers can be easily created by following the steps below:

Please replace the letters in **green** with the part no. code as per your requirements.

1. Components' grounding:

S	=	not grounded
A	=	grounded

2. Choose the arm version required:

910G / 910GS	=	flush mount
224G / 230G / 224GS / 230GS	=	support arm pivot
33030G / 33030GS	=	extended pivot arm light duty
H33030G / H33030GS	=	extended pivot arm heavy duty
431G / 431GS	=	height adjustable arm
73130G / 73130GS	=	extended height adjustable arm
531G / 531GS	=	high position height adjustable arm
83130G / 83130GS	=	extended high position height adjustable arm

3. Specify where you want to mount the arm and select the fixing (Y):

1V	=	pole ø 23 - 40 mm / 0.9 - 1.6"
1	=	pole ø 25 mm / 1.0"
1B	=	pole ø 30 mm / 1.2"
1F	=	pole ø 35 mm / 1.4"
1C	=	pole ø 38 mm / 1.5"
2F	=	horizontal rails 10 x 25, 8 x 35, 10 x 30, 10 x 50 mm
3	=	vertical rail 10 x 25 mm
4S	=	vertical wall channel
9	=	Löwenstein Medical / ITD / mounting channel
13	=	Mindray medical supply units
17	=	Amico pendants & headwalls

How to create part numbers

4. Choose the colour requested for decorative parts (C):

gr	=	RAL 7024 (graphite grey)
bl	=	RAL 5013 (cobalt blue)

5. Choose the colour requested for the aluminium die-cast parts (D):

1	=	RAL 9016 (traffic white)
2	=	RAL 9002 (grey white)

6. Specify the spring depending on the weight of your monitor/ device - only with height adjustable arms (X):

60	=	2 - 5 kg / 4.4 - 11 lbs
95	=	4 - 9 kg / 8.8 - 19.8 lbs
110	=	4.5 - 12.5 kg / 9.9 - 27.6 lbs
130	=	7.5 - 13 kg / 16.5 - 28.7 lbs
160	=	12.5 - 19 kg / 27.6 - 41.9 lbs
175	=	13 - 16 kg / 28.7 - 35.3 lbs
220	=	16 - 22 kg / 35.3 - 48.5 lbs

Example

01 Grounding		02 Arm version		03 Fixing		04 & 05 Colour		06 Spring tension
A	-	431G	-	Y	.	CD	-	X
A	-	431G	-	4S	.	gr1	-	95

Product highlights

Non-height adjustable arm

Pivot and extended pivot arms offer **full cable integration** and easy cleaning for **infection control** in all critical care areas in the clinical environment. Our **5-years warranty** guarantees reliable, trouble-free and maintenance-free operation. The swivel function at the fixing point as well as the **tilt/ swivel function** of the monitor adaptation offers optimum positioning of your device. The extended pivot arm is a two-part arm with central pivot point between the two arm sections. All bottom-mounted monitors have the safety feature of a **20° tilt stop**.

Height adjustable arm


The height adjustable arm offers **full cable integration**, one-wipe clean surfaces and our **5-year warranty** on the entirety of the product. Vertically and horizontally adjustable, the arm offers optimal positioning of your device. A constant viewing angle when vertically adjusted is guaranteed by a parallel linkage. **Weight-optimized springs counter-balance** the monitor without making any adjustments. All bottom-mounted monitors have the safety feature of a **20° tilt stop**. Individual rotational stops may be integrated in any joint.

To provide the perfect combination of arm and monitor, the CIM articulated height adjustable arm has a simple quarter-turn lock to hold the arm safely in position when removing the monitor. Furthermore the arm is equipped with the unique safety feature of a slow-release gas spring. Should anybody try to remove the monitor before locking the arm in position, the arm will only rise at a slow and even tempo.

SAFETY

Flush mount

MD

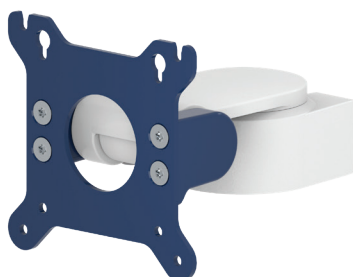

Monitor adaption

VESA 75/100

Colours to choose CD

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

18 kg / 39.7 lbs

Max. load monitor adaptation

18 kg / 39.7 lbs

Product weight

1 kg / 2.2 lbs VESA 75/100
1.1 kg / 2.4 lbs VESA 75/100 rotatable

Part number

S-910G-Y.CD VESA 75/100
S-910GS-Y.CD VESA 75/100 rotatable

Fixings Y

details page 31




4S

Adapter for vertical wall channel



2F

Horizontal rail clamp

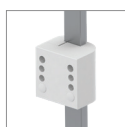

3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp



13

Adapter for Mindray medical supply units




17

Adapter for Amico pendants & headwalls


Flush mount

Technical data

SIDE VIEW

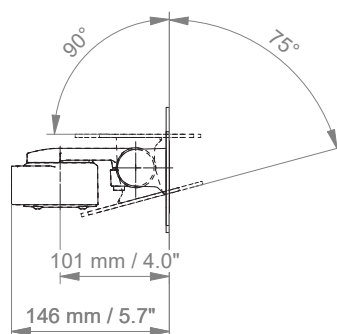


TOP VIEW

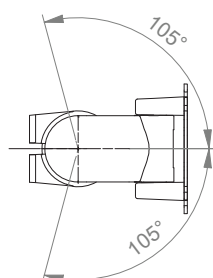

Support arm pivot

MD

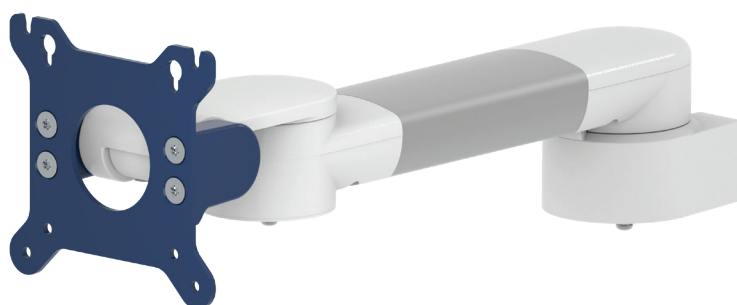

Monitor adaption

VESA 75/100

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

18 kg / 39.7 lbs

Max. load monitor adaptation

18 kg / 39.7 lbs

Product weight

arm type 240	arm type 300	
1.8 kg / 4.0 lbs	1.9 kg / 4.2 lbs	VESA 75/100
1.9 kg / 4.2 lbs	2.1 kg / 4.6 lbs	VESA 75/100 rotatable

Part number

arm type 240	arm type 300	
A-224G- Y.CD	A-230G- Y.CD	VESA 75/100
A-224GS- Y.CD	A-230GS- Y.CD	VESA 75/100 rotatable

Fixings **Y**

details page 31



4S

Adapter for vertical wall channel



2F

Horizontal rail clamp



3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp




13

Adapter for Mindray medical supply units



17

Adapter for Amico pendants & headwalls

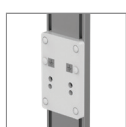


9

Adapter for Löwenstein Medical / ITD / mounting channel


Support arm pivot

Technical data - arm type 300

SIDE VIEW

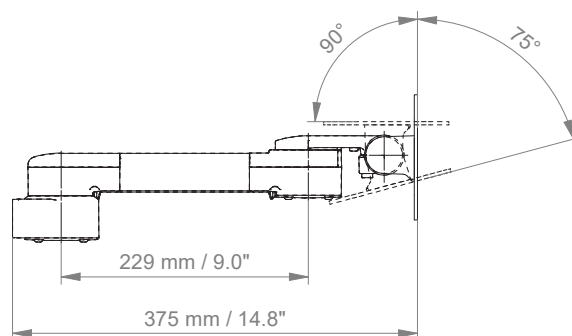


TOP VIEW

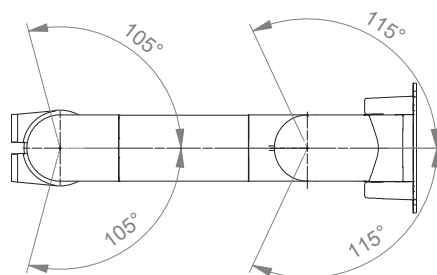

Extended pivot arm light duty

MD

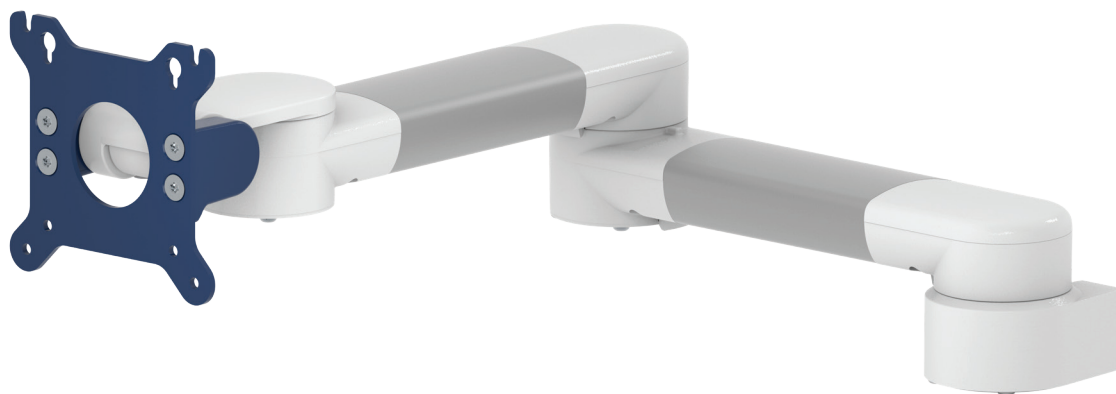

Monitor adaption

VESA 75/100

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

8 kg / 17.6 lbs

Max. load monitor adaptation

8 kg / 17.6 lbs

Product weight

2.9 kg / 6.4 lbs VESA 75/100
3.0 kg / 6.6 lbs VESA 75/100 rotatable

Part number

A-33030G-**Y.CD** VESA 75/100
A-33030GS-**Y.CD** VESA 75/100 rotatable

Fixings **Y**

details page 31




4S

Adapter for
vertical wall channel



2F

Horizontal rail clamp

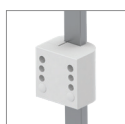

3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp




13

Adapter for Mindray
medical supply units



17

Adapter for Amico
pendants & headwalls

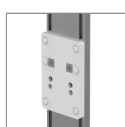


9

Adapter for Löwenstein Medical /
ITD / mounting channel


Extended pivot arm light duty

Technical data

SIDE VIEW

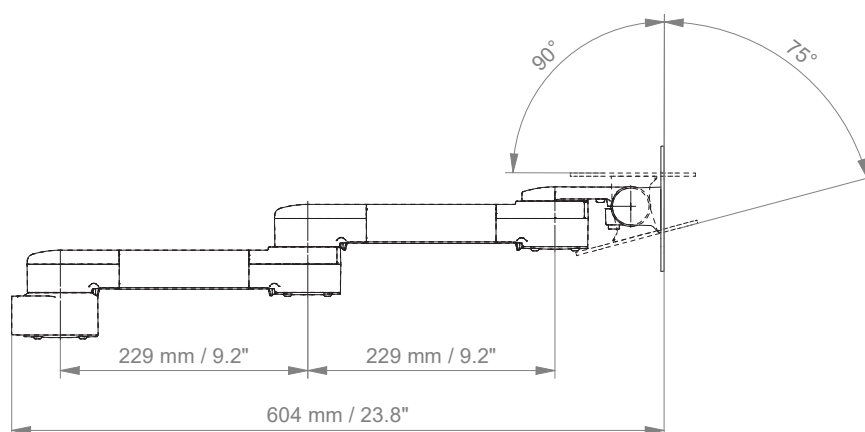

TOP VIEW

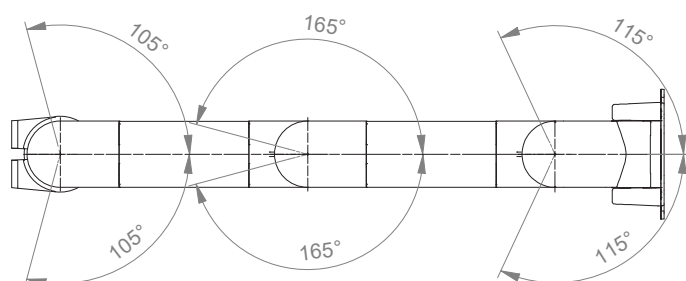

Extended pivot arm heavy duty

MD



Monitor adaption

VESA 75/100

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

18 kg / 39.7 lbs

Max. load monitor adaptation

18 kg / 39.7 lbs

Product weight

3.1 kg / 6.8 lbs VESA 75/100
3.2 kg / 7.1 lbs VESA 75/100 rotatable

Part number

A-H33030G-**Y.CD** VESA 75/100
A-H33030GS-**Y.CD** VESA 75/100 rotatable

Fixings **Y**

details page 31




4S

Adapter for
vertical wall channel



2F

Horizontal rail clamp

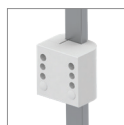

3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp




13

Adapter for Mindray
medical supply units



17

Adapter for Amico
pendants & headwalls

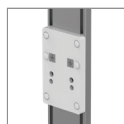


9

Adapter for Löwenstein Medical /
ITD / mounting channel


Extended pivot arm heavy duty

Technical data


SIDE VIEW

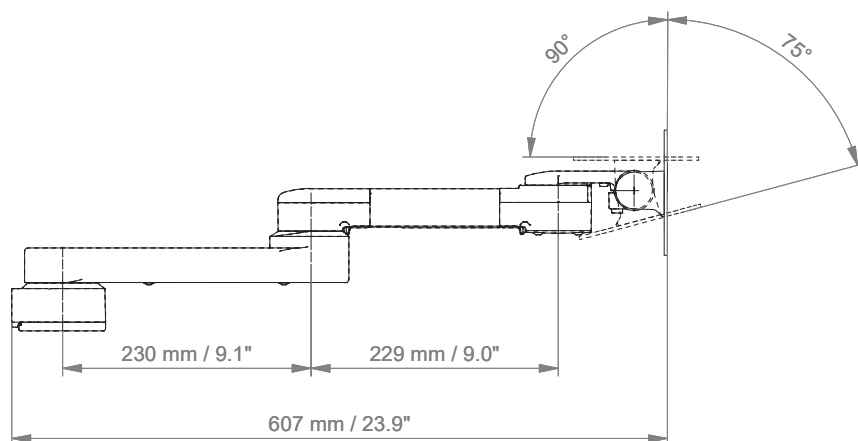

TOP VIEW

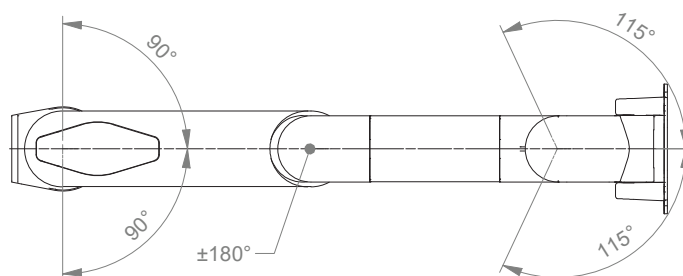

Height adjustable arm

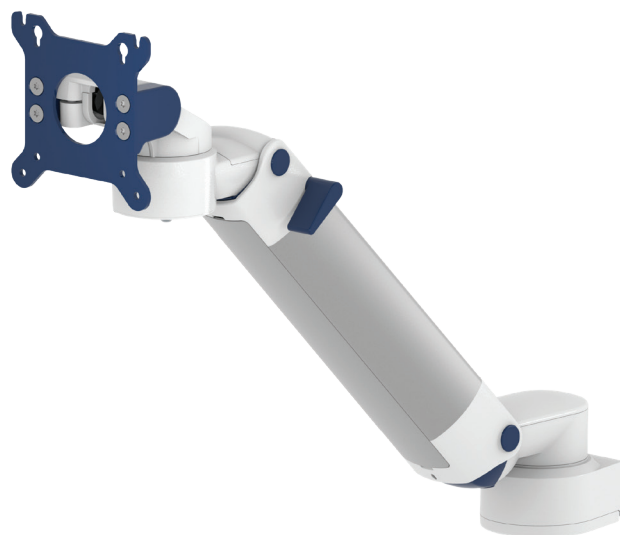MD

Monitor adaption

VESA 75/100

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Colours to choose **CD**

decorative parts	RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts	RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

22 kg / 48.5 lbs

Max. load monitor adaption

18 kg / 39.7 lbs

Product weight

3.8 kg / 8.4 lbs	VESA 75/100
3.9 kg / 8.6 lbs	VESA 75/100 rotatable

Balancing range **X**

600 N (2-5 kg / 4.4-11 lbs); 950 N (4-9 kg / 8.8-19.8 lbs); 1100 N (4.5-12.5 kg / 9.9-27.6 lbs); 1300 N (7.5-13 kg / 16.5-28.7 lbs); 1600 N (12.5-19 kg / 27.6-41.9 lbs); 1750 N (13-16 kg / 28.7-35.3 lbs); 2200 N (16-22 kg / 35.3-48.5 lbs)

Part number

A-431G-Y.CD_X	VESA 75/100
A-431GS-Y.CD_X	VESA 75/100 rotatable

Fixings **Y**


details page 31


4S

Adapter for vertical wall channel


2F

Horizontal rail clamp

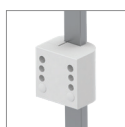
3

Vertical rail clamp


1V

Universal pole clamp


1
1B
1F
1C


Fixed pole clamp


13

Adapter for Mindray medical supply units


17

Adapter for Amico pendants & headwalls

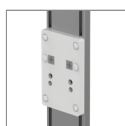

9

Adapter for Löwenstein Medical / ITD / mounting channel


Height adjustable arm

Technical data

SIDE VIEW

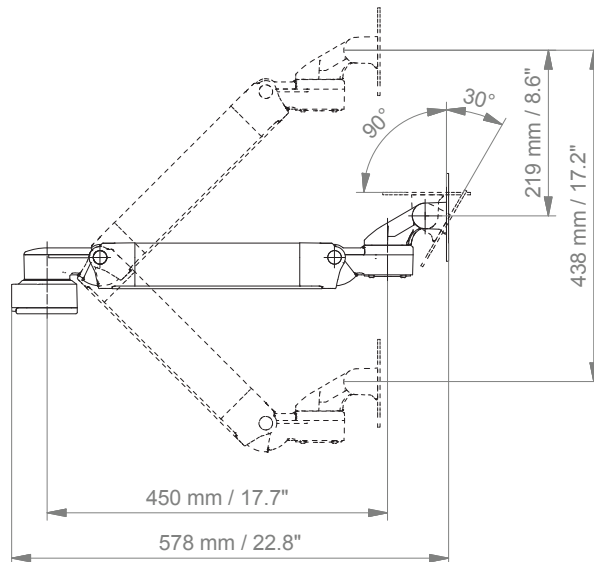

TOP VIEW

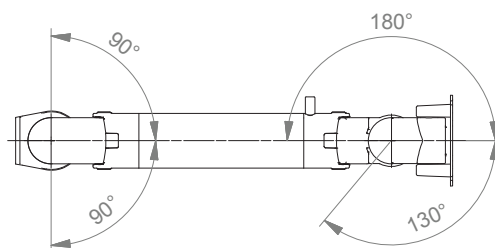

Extended height adjustable arm

MD



Monitor adaption

VESA 75/100

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

22 kg / 48.5 lbs

Max. load monitor adaption

18 kg / 39.7 lbs

Product weight

4.8 kg / 10.6 lbs VESA 75/100
4.9 kg / 10.8 lbs VESA 75/100 rotatable

Balancing range **X**

600 N (2-5 kg / 4.4-11 lbs); 950 N (4-9 kg / 8.8-19.8 lbs); 1100 N (4.5-12.5 kg / 9.9-27.6 lbs); 1300 N (7.5-13 kg / 16.5-28.7 lbs); 1600 N (12.5-19 kg / 27.6-41.9 lbs); 1750 N (13-16 kg / 28.7-35.3 lbs); 2200 N (16-22 kg / 35.3-48.5 lbs)

Part number

A-73130G-**Y.CD_X** VESA 75/100
A-73130GS-**Y.CD_X** VESA 75/100 rotatable

Fixings **Y**

details page 31




4S

Adapter for vertical wall channel



2F

Horizontal rail clamp

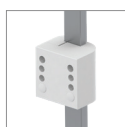

3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp




13

Adapter for Mindray medical supply units



17

Adapter for Amico pendants & headwalls

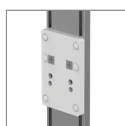


9

Adapter for Löwenstein Medical / ITD / mounting channel


Extended height adjustable arm

Technical data

SIDE VIEW

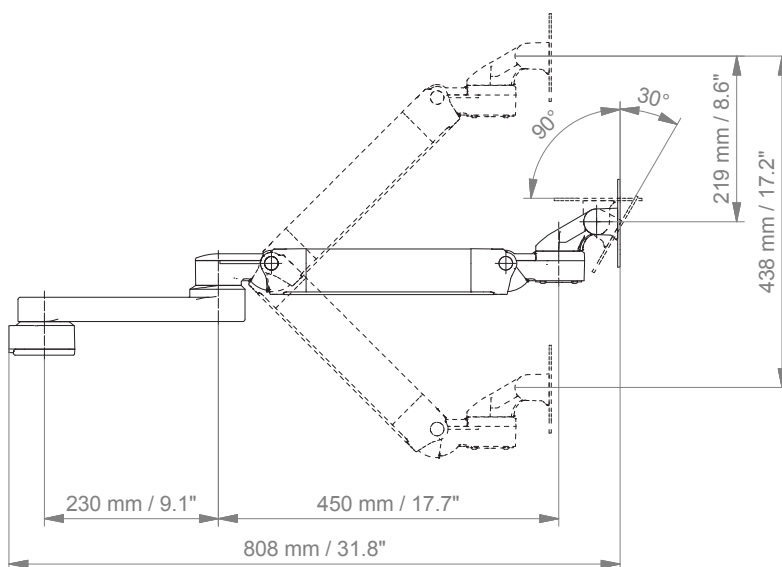


TOP VIEW

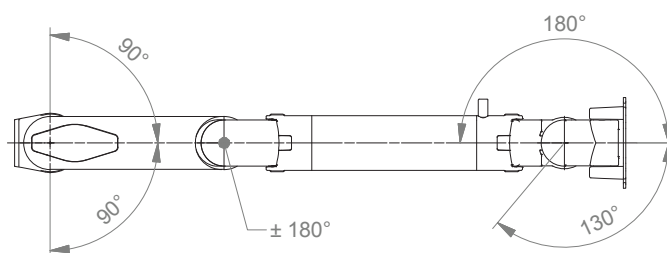

High position height adjustable arm

MD

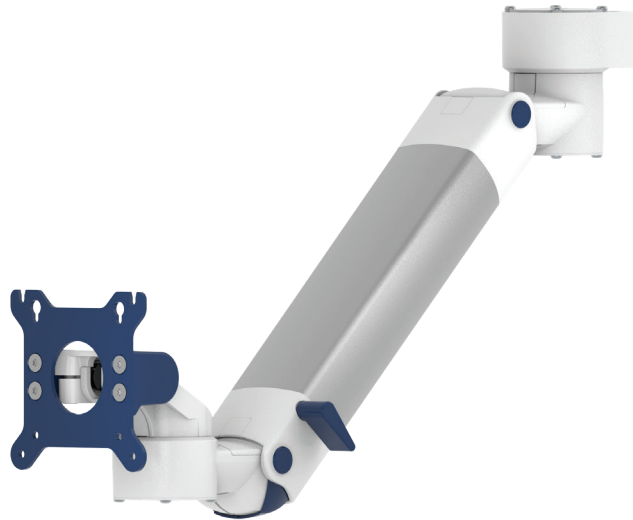

Monitor adaption

VESA 75/100

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

22 kg / 48.5 lbs

Max. load monitor adaptation

18 kg / 39.7 lbs

Product weight

4.1 kg / 9.0 lbs VESA 75/100
4.3 kg / 9.5 lbs VESA 75/100 rotatable

Balancing range **X**

600 N (2-5 kg / 4.4-11 lbs); 950 N (4-9 kg / 8.8-19.8 lbs); 1100 N (4.5-12.5 kg / 9.9-27.6 lbs); 1300 N (7.5-13 kg / 16.5-28.7 lbs); 1600 N (12.5-19 kg / 27.6-41.9 lbs); 1750 N (13-16 kg / 28.7-35.3 lbs); 2200 N (16-22 kg / 35.3-48.5 lbs)

Part number

A-531G-**Y**.**CD**_X VESA 75/100
A-531GS-**Y**.**CD**_X VESA 75/100 rotatable

Fixings **Y**

details page 31



4S

Adapter for vertical wall channel



3

Vertical rail clamp



1V

Universal pole clamp



1
1B
1F
1C

Fixed pole clamp



13

Adapter for Mindray medical supply units




17

Adapter for Amico pendants & headwalls


High position height adjustable arm

Technical data

SIDE VIEW

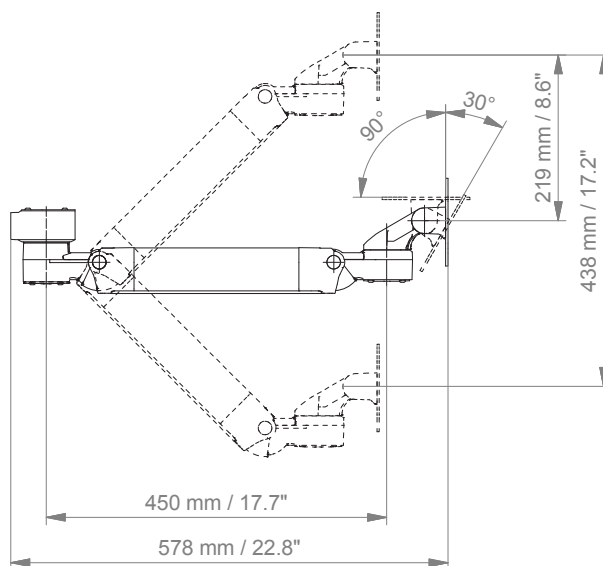


TOP VIEW

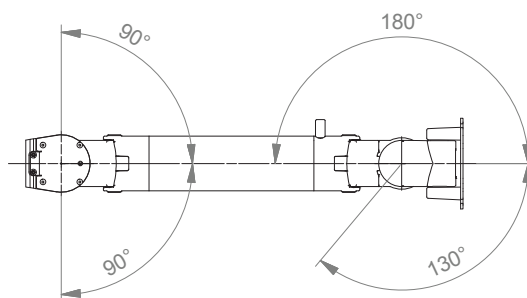

Extended high position height adjustable arm

MD

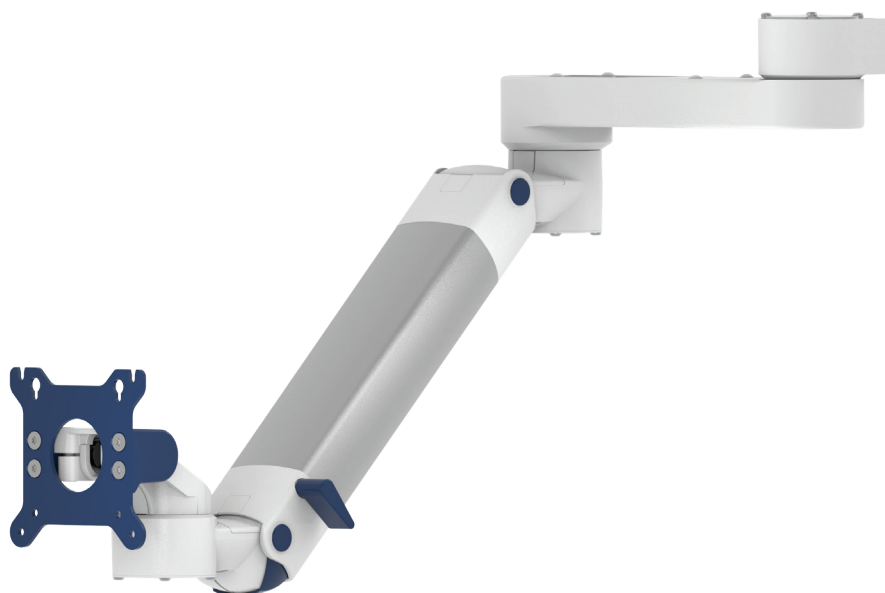

Monitor adaption

VESA 75/100

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load capacity support arm

22 kg / 48.5 lbs

Max. load monitor adaptation

18 kg / 39.7 lbs

Product weight

5.3 kg / 11.7 lbs

VESA 75/100

5.4 kg / 11.9 lbs

VESA 75/100 rotatable

Balancing range **X**

600 N (2-5 kg / 4.4-11 lbs); 950 N (4-9 kg / 8.8-19.8 lbs); 1100 N (4.5-12.5 kg / 9.9-27.6 lbs); 1300 N (7.5-13 kg / 16.5-28.7 lbs); 1600 N (12.5-19 kg / 27.6-41.9 lbs); 1750 N (13-16 kg / 28.7-35.3 lbs); 2200 N (16-22 kg / 35.3-48.5 lbs)

Part number

A-83130G-**Y**.**CD**_X VESA 75/100

A-83130GS-**Y**.**CD**_X VESA 75/100 rotatable

Fixings **Y**

details page 31



4S

Adapter for vertical wall channel



3

Vertical rail clamp



1V

Universal pole clamp



1

1B

1F

1C

Fixed pole clamp



13

Adapter for Mindray medical supply units




17

Adapter for Amico pendants & headwalls


Extended high position height adjustable arm

Technical data

SIDE VIEW

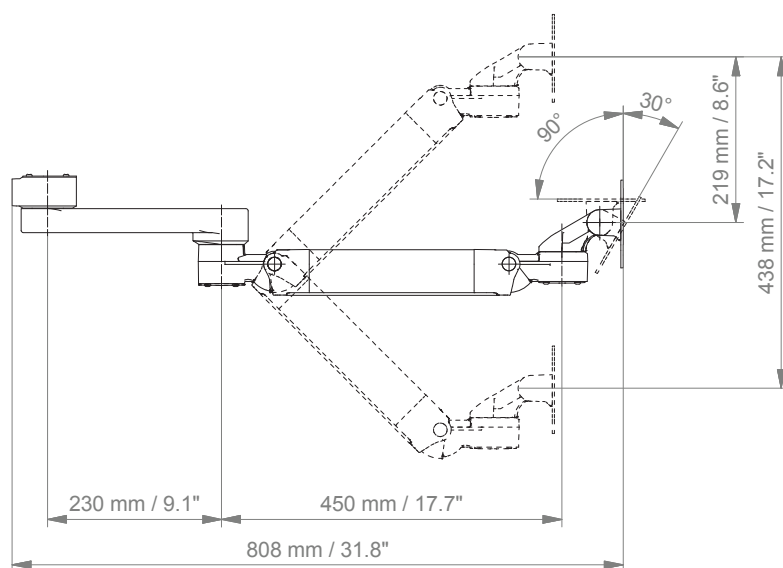

TOP VIEW

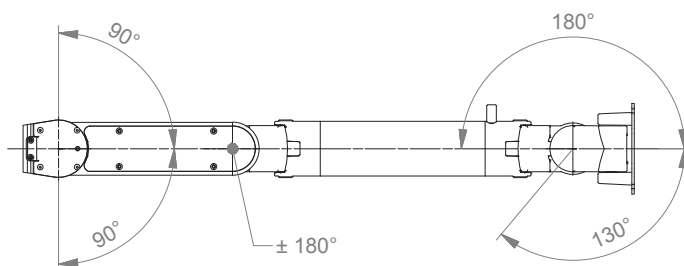

Counter top

MD



Monitor adaption

VESA 75/100

Colours to choose CD

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load monitor adaptation

18 kg / 39.7 lbs

Part number

Counter top with base plate to mount to desk or even surfaces

T-HG-60.CD_220 height 220 mm / 8.7"
T-HG-60.CD_280 height 280 mm / 11"
T-HG-60.CD_340 height 340 mm / 13.4"

Counter top free-standing (see illustration)

T-HG-62.CD_220 height 220 mm / 8.7"
T-HG-62.CD_280 height 280 mm / 11"
T-HG-62.CD_340 height 340 mm / 13.4"

Counter top without base plate

T-HG-69.CD_220 height 220 mm / 8.7"
T-HG-69.CD_280 height 280 mm / 11"
T-HG-69.CD_340 height 340 mm / 13.4"




All free-standing counter tops
have to be tip-tested for stability
depending on the weight of the
monitor!


Counter top

Technical data - free standing

SIDE VIEW

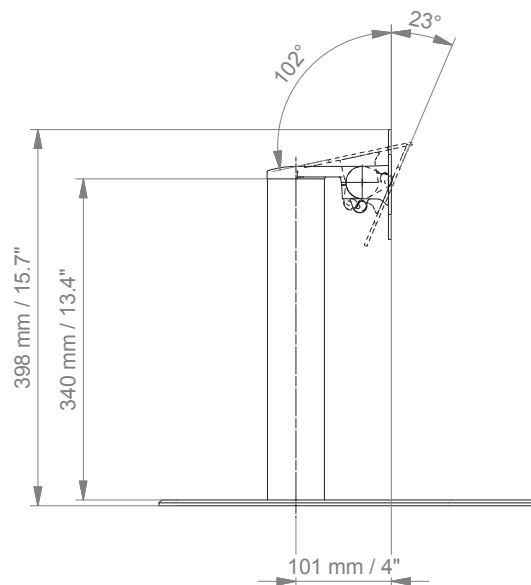

TOP VIEW

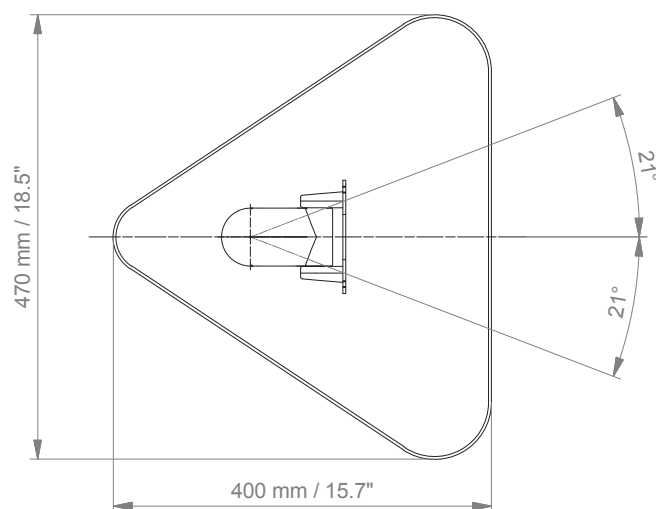

Counter top for deep monitors

MD



Monitor adaption

CL-bracket, VESA 75/100

Colours to choose **CD**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white / RAL 9002 grey white

Max. load monitor adaptation

18 kg / 39.7 lbs

Part number

Counter top with base plate
to mount to desk or even surfaces (see illustration)
T-HK-60.**CD**_50

Counter top free-standing
T-HK-62.**CD**_50




All free-standing counter tops
have to be tip-tested for stability
depending on the weight of the
monitor!


Counter top for deep monitors

Technical data - with base plate

SIDE VIEW

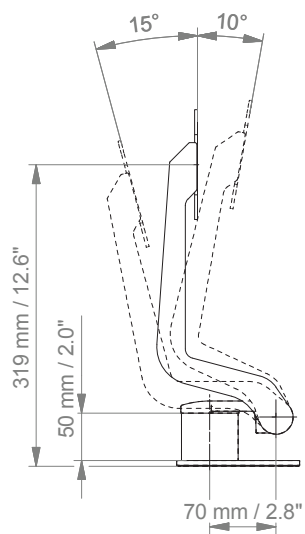

TOP VIEW

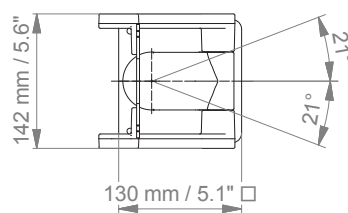

Cart light duty

MD

Example configuration



Monitor adaption

CL-bracket, VESA 75/100

Accessories

5-star base aluminum base with counter weight
double wheel castors, ø 100 mm / 3.9", 2 lockable
column to accept components and accessories on opposite sides

Colours to choose **CD**

decorative parts RAL 5013 blue / RAL 7024 graphite grey
base RAL 9016 pure white
handle RAL 9016 pure white / RAL 5013 blue / RAL 7024 graphite grey

Max. load capacity

10 kg / 22 lbs

Optional accessories

basket big 300 x 200 x 180 mm / 11.8 x 7.9 x 7.1"
basket small 200 x 150 x 100 mm / 7.9 x 5.9 x 3.9"
handle


Part number

cart	R-002K. CD
basket big	Z-51-4Z.1
basket small	Z-52-4Z.1
handle	Z-HD-4. C1


Cart light duty

Technical data

SIDE VIEW

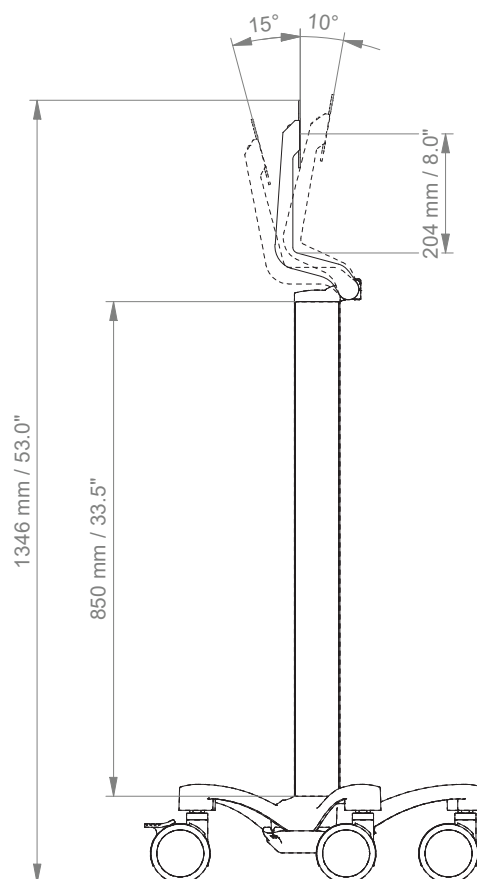

TOP VIEW

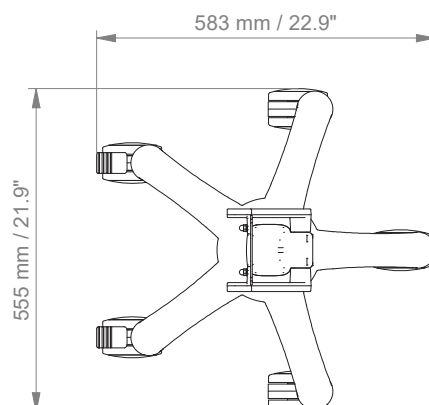

Cart light duty kit

MD

Kit configuration



Monitor adaption

CL-bracket, VESA 75/100

Cart light duty

5-star base aluminum base with counter weight
double wheel castors, ø 100 mm / 3.9", 2 lockable
column to accept components and accessories on opposite sides

Colours

decorative parts RAL 5013 cobalt blue
base RAL 9016 traffic white
handle RAL 7024 graphite grey

Max. load monitor adaptation

10 kg / 22 lbs

Max. load keyboard tray

1 kg / 2.2 lbs

Kit number

R-002K.gr1_286490


Consisting of the following parts:

R-002K.gr1	cart light duty
Z-HD-4.gr1	handle for carts
Z-EF-4V.gr1_50S	keyboard tray magnetic with fold-up function 500 x 200 mm/ 19.7 x 7.9"


Cart light duty kit

Technical data

SIDE VIEW

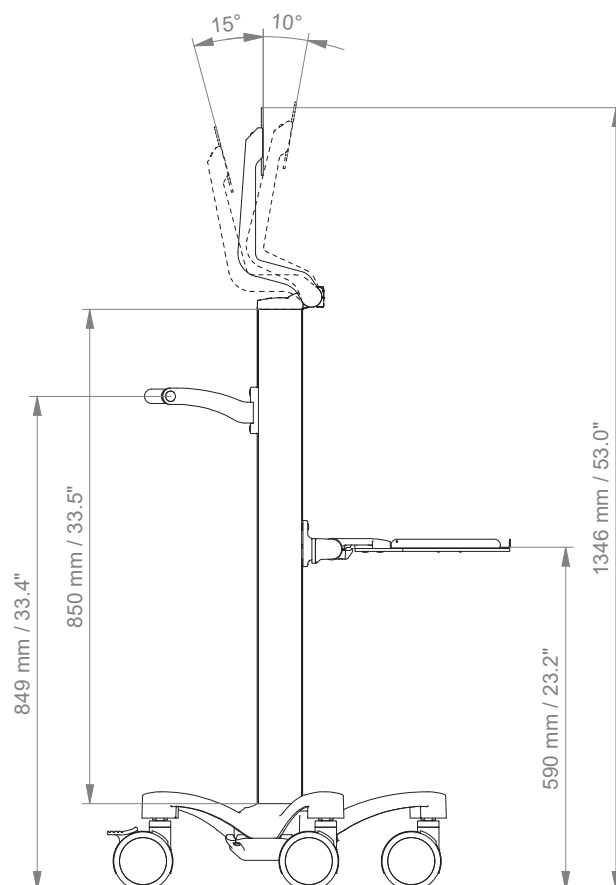

TOP VIEW

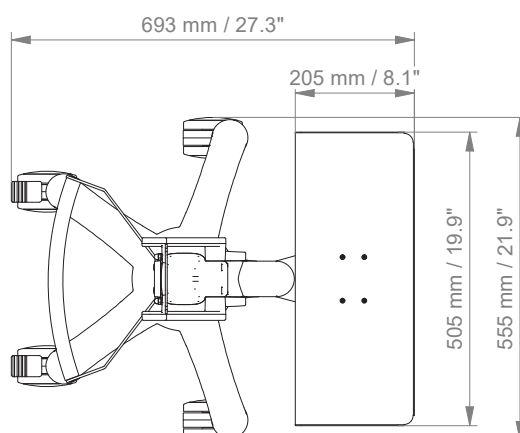

Fixings Y



Illustration	Description	Notes	Weight	Part no.
	Adapter for vertical wall channel	To mount to CIM wall channel (compatible with GCX channel) Please order separately: CIM wall channel, see catalogue accessoires	0.6 kg / 1.3 lbs	4S
	Horizontal rail clamp	To mount to horizontal European rails 10 x 25, 8 x 35, 10 x 30, 10 x 50 mm	1.1 kg / 2.4 lbs	2F
	Vertical rail clamp	To mount to vertical European rail 10 x 25 mm	0.5 kg / 1.1 lbs	3
	Universal pole clamp	V-block to mount to pole $\varnothing$ 23 - 40 mm / 0.9 - 1.6"	0.6 kg / 1.3 lbs	1V
	Fixed pole clamp	To mount to pole $\varnothing$ 25 mm / 1.0" To mount to pole $\varnothing$ 30 mm / 1.2" To mount to pole $\varnothing$ 35 mm / 1.4" To mount to pole $\varnothing$ 38 mm / 1.5"	0.3 kg / 0.7 lbs 0.4 kg / 0.9 lbs 0.4 kg / 0.9 lbs 0.4 kg / 0.9 lbs	1 1B 1F 1C
	Adapter medical supply units	For Mindray medical supply units	0.8 kg / 1.8 lbs	13
	Adapter pendants & headwalls	For Amico pendants & headwalls	0.6 kg / 1.3 lbs	17
	Adapter pendants & headwalls	For Löwenstein Medical / ITD / mounting channels	0.4 kg / 0.9 lbs	9



CIM med GmbH

Euro-Industriepark
Margot-Kalinke-Str. 9
80939 Munich / Germany

T +49 89. 978 94 08-0

CIM med Inc.

2024 E Independence Commerce Dr
Matthews, NC 28105 / USA

T +1 855 895 4722

info@cim-med.com
www.cim-med.com

